

(19)

Europäisches Patentamt
European Patent Office
Office européen des brevets



AC

(11)

EP 1 439 191 A1

(12)

EUROPEAN PATENT APPLICATION

(43) Date of publication:
21.07.2004 Bulletin 2004/30

(51) Int Cl.7: **C07K 14/52, C12N 15/19**

(21) Application number: **03100064.9**

(22) Date of filing: **15.01.2003**

(84) Designated Contracting States:
AT BE BG CH CY CZ DE DK EE ES FI FR GB GR
HU IE IT LI LU MC NL PT SE SI SK TR
Designated Extension States:
AL LT LV MK RO

(72) Inventors:
• **Capponi, Luciano**
00050 Fiumicino (IT)
• **Rossi, Mara**
00179 Roma (IT)

(71) Applicant: **ARES TRADING S.A.**
CH-2028 Vaumarcus (CH)

(74) Representative: **Serono International S.A.**
12, Chemin des Aulx
1228 Plan-les-Ouates (CH)

(54) **Process for the purification of bacterially expressed proteins**

(57) This patent application relates to a process for the purification of proteins expressed in prokaryotic cells, which comprises interposing a reverse phase chromatography step between the extraction/denaturation and the renaturation steps.

EP 1 439 191 A1

DescriptionFIELD OF THE INVENTION

5 **[0001]** This invention relates to a process for the purification of proteins expressed in prokaryotic cells

BACKGROUND OF THE INVENTION

10 **[0002]** Recombinant DNA technology makes it possible to produce large amounts of desired proteins. Bacteria are a particularly convenient source for recombinant protein production for purification purposes. They can be grown in very large quantities under well-defined conditions, and they are relatively easy to break open for extraction. Most importantly, molecular cloning techniques allow high levels of expression in bacteria of almost any protein from any organism, bacterial or otherwise. Unfortunately, bacterially expressed proteins are often difficult to purify because of their tendency to precipitate within the cell. The precipitated protein forms inclusion bodies: dense, granular structures that are distributed throughout the cytoplasm. Inclusion body formation is especially common for nonbacterial proteins, although even native bacterial proteins show a tendency to aggregate when they are expressed at very high levels. It is not known exactly how they are formed, but it is thought that the protein is partially or incorrectly folded. The major disadvantage of inclusion bodies is that extraction of the protein of interest generally requires the use of denaturing agents. This can cause problems when a native folded protein is required, because refolding methods are not always 100 percent effective and can be difficult to scale-up.

20 **[0003]** The advantage of inclusion bodies is that they generally allow greater levels of expression, and they can be separated easily from a large proportion of bacterial cytoplasmic proteins by centrifugation, providing an effective purification step.

25 **[0004]** For bacterial E. coli lysis, techniques such as mechanical lysis, sonication, enzymatic lysis and detergent lysis are available. These techniques, however, either require special instrumentation or have certain limitations. Most importantly, these methods can neither fully recover soluble recombinant proteins, nor regain inclusion bodies. Thus, the yield of recombinant proteins can be very low. It is necessary to explore alternative methods in order to improve protein yields by recovering both the soluble fraction and insoluble fraction of the recombinant protein.

30 DESCRIPTION OF THE INVENTION

35 **[0005]** According to the present invention, a Reverse Phase Chromatography (RPC) step is interposed between the extraction/denaturation step and the subsequent refolding step, resulting in an increased yield of the protein and other advantages that will become apparent in the following description. It is believed that the removal of most cell-derived impurities in the RPC step helps in obtaining a higher yield in the subsequent refolding step, although the invention should be considered independently from this hypothesis. Compounds stick to reverse phase HPLC columns in high aqueous mobile phase and are eluted from RP HPLC columns with high organic mobile phase. In RP HPLC compounds are separated based on their hydrophobic character. Since columns are tubular, column dimensions usually take the following format, internal diameter X length (for example 4.6mm X 250mm).

40 **[0006]** The stationary phase is generally made up of hydrophobic alkyl chains ($-\text{CH}_2-\text{CH}_2-\text{CH}_2-\text{CH}_3$) that interact with the analyte.

45 **[0007]** In fact, media for RPC are typically highly substituted with hydrophobic ligands and the binding of substances to RPC media is usually very strong and requires organic solvents for elution.

50 **[0008]** SOURCETM (Amersham) RPC can be used at full pH range. SOURCETM RPC are designed for fast, high performance preparative separations of bio-molecules such as proteins, peptides and oligonucleotides. The media have matrices based on rigid, polystyrene/divinyl benzene, with monosized beads of diameters 15 μm or 30 μm respectively.

55 **[0009]** Pore size distribution is controlled and reproducible. The wide pH stability and high capacity make SOURCETM RPC media an interesting alternative to silica-based media. The high chemical stability of the matrix offers unmatched flexibility in choice of running and cleaning conditions.

60 **[0010]** The reverse phase solvents are by convention installed on the HPLC channels A and B. The A solvent by convention is the aqueous solvent (water) and the B solvent by convention is the organic solvent (acetonitrile, methanol, propanol). It is important to follow this convention since the terms A and B are commonly used to refer to the aqueous and organic solvents respectively. The A solvent is generally HPLC grade water with 0.1% acid. The B solvent is generally an HPLC grade organic solvent such as acetonitrile or methanol with 0.1% acid. The acid is used to improve the chromatographic peak shape and to provide a source of protons in reverse phase LC/MS. The acids most commonly used are formic acid, trifluoroacetic acid, and acetic acid.

65 **[0011]** As already said, the process of the invention represents an improvement in this process consisting in that

a Reverse Phase Chromatography (RPC) step is interposed between the extraction/denaturation step and the subsequent refolding step. However, all the steps normally carried out in a classical process for purifying a protein recombinantly produced in prokaryotic host cells are reported here below just for completeness.

[0012] Normally a classical process for purifying a protein produced in prokaryotic cells also includes the following steps:

- Preparation of cell lysates from the host cells
- Isolation of inclusion bodies
- Solubilization of the aggregated proteins in the inclusion bodies/denaturation step
- Refolding/renaturation of the solubilized proteins.

Preparation of cell lysates from the host cells

[0013] The protein expressed in prokaryotic cells needs first to be extracted from the host cells in which it has been expressed. A variety of methods are available to lyse cells. Which of these methods should be used in a specific case depends on the type of host cells and on the quantity of cells to be lysed.

[0014] The first choice to make is that of the nature and the pH of the buffer system we want to use. This depends on:

- the stability of the target protein with respect to pH and the buffering compound.
- the purification procedure.

[0015] To avoid time and protein loss caused by an additional buffer exchange step, it is advisable to choose a buffer that is compatible with the first chromatography step (see chromatography). Buffers and their pH ranges are listed in Table 1. The most used buffers are phosphate, Tris - HCl and HEPES NaOH. They are normally used at concentrations of 20-50 mM.

Table 1

Buffer	pH range
Citric acid - NaOH	2.2 - 6.5
Sodium citrate - citric acid	3.0 - 6.2
Sodium acetate - acetic acid	3.6 - 5.6
Cacodylic acid sodium salt - HCl	5.0 - 7.4
MES-NaOH	5.6 - 6.8
Sodium dihydrogen phosphate - disodium hydrogen phosphate	5.8 - 8.0
Imidazole - HCl	6.2 - 7.8
MOPS - KOH	6.6 - 7.8
Triethanolamine hydrochloride - NaOH	6.8 - 8.8
Tris - HCl	7.0 - 9.0
HEPES - NaOH	7.2 - 8.2
Tricine - NaOH	7.6 - 8.6
Sodium tetraborate - boric acid	7.6 - 9.2
Bicine - NaOH	7.7 - 8.9
Glycine - NaOH	8.6 - 10.6

[0016] Depending on the target protein, it may be necessary to add compounds to the lysis buffer to improve the stability of the target protein and to keep the protein in solution. The most used additives, their effective concentrations, and their general purpose are listed in Table 2.

Table 2

Class of Additive	Example	Concentration	Purpose
Salts	NaCl, KCl, $(\text{NH}_4)_2\text{SO}_4$	50-150 mM	maintain ionic strength of medium
Detergents	Deoxycholate, Triton X-100	0.1-1%	solubilization of poorly soluble proteins

Table 2 (continued)

Class of Additive	Example	Concentration	Purpose
Glycerol		5-10%	stabilization
Glucose or sucrose		25 mM	Stabilize lysosomal membranes, reduce protease release
Metal chelators	EDTA, EGTA	1 mM	reduce oxidation damage, chelate metal ions
Reducing agents	2-Mercaptoethanol, DTT	0.05%	reduce oxidation
Ligands, metal ions	Mg ²⁺ , ATP, GTP	1-10 mM	stabilization

[0017] An important class of additives is the protease inhibitors. In general cell disruption leads to the release of proteolytic enzymes, which could lower the overall yield. To control this undesirable proteolysis it may be necessary to add a cocktail of protease inhibitors to the cell suspension. Since many of these compounds are not very stable in aqueous solutions it is important that they are added to the lysis buffer from a stock solution in an organic solvent (methanol, ethanol, isopropanol, or DMSO) immediately before use.

[0018] Although it is impossible to give a general lysis buffer, a good starting buffer would be:

50 mM Tris-HCl pH 7.5
100 mM NaCl
1 mM DTT (Dithiothreitol)
5% glycerol (possibly).

[0019] Most lysis methods cause the release of nucleic acids (DNA and RNA). These have to be removed because they can cause viscosity problems or due to their interference with subsequent chromatographic steps. Different methods exist:

- Enzymatic digestion by the addition of DNase I (1 µg/ml) to the cell lysate. The mixture is incubated on ice for 10-15 min.
- Mechanical breakdown by shearing during sonication. When the French Pressure Cell is used it is advisable to add DNase to the cell suspension.
- Precipitation by treatment with polyethyleneimine (0.1% (w/v)) or protamine sulphate (1% (w/v)) followed by centrifugation. Add the precipitants to the cell lysate and incubate the solution for 30 min at 4°C.

[0020] Different methods are used for the preparation of cell lysates from *E. coli* cells.

[0021] Sonication is the most popular technique for lysing small quantities of cells (1-6 L of cell culture). Cells are lysed by liquid shear and cavitation. DNA is also sheared during sonication, so it is not necessary to add DNase to the cell suspension. The main problem is controlling the temperature. This is addressed by keeping the suspension on ice and using a number of short pulses (5-10 sec) with pauses (10-30 sec) to reestablish a low temperature. For cell quantities larger than 50 g the method is of limited value because of the difficulty in maintaining low temperatures and the long sonication times needed to reach adequate lysis.

[0022] Homogenizers are the most common devices to lyse bacteria. The presses lyse cells by pressurizing the cell suspension and suddenly releasing the pressure. This creates a liquid shear capable of lysing cells. Typical operating pressures for the older type of homogenizers, the French press and Manton -Gaulin homogenizer, are 6000-10,000 psi. Multiple (2-3) passes are generally required to achieve adequate lysis. The high operating pressures, however, result in a rise in operating temperatures.

[0023] Therefore, pressure cells are cooled (4°C) prior to use. In addition to temperature control, care should be taken to avoid inactivating proteins by foaming.

[0024] Modern homogenizers are often continuous and can be operated at higher pressures. At EMBL we have an Avestin Emulsiflex-C5 that efficiently lyses *E. coli* cells in one passage at 15,000 psi (100 MPa).

[0025] Enzymatic lysis is based on the digestion of the peptidoglycan layer of the bacterial cell wall by lysozyme. Gram-negative bacteria, however, have an outer membrane that is external to the cell wall and needs to be permeabilized to expose the peptidoglycan layer. Tris, often used as a buffer in lysis methods, effectively permeabilizes outer membranes. This effect can be enhanced by the addition of EDTA (1 mM). EDTA chelates the magnesium ions that stabilize membranes. During cell lysis often a lot of DNA is liberated and it becomes necessary to add DNase (1 mg/

ml) to reduce the viscosity of the preparation. Enzymatic cell lysis can be carried out on any scale but for large-scale preparations the lysosyme and DNase can get expensive. To increase the level of cell lysis the solution can be sonicated (see above).

[0026] An alternative lysis method is to freeze the cells directly in liquid nitrogen and ground the frozen cells to a powder using a mortar and pestle that are chilled with liquid nitrogen. The powder can be stored indefinitely at -80°C and the cell lysate can be prepared by adding the powder to 5 volumes of buffer.

[0027] The extent of the solubilization and the stability of the solubilized protein depend on the detergent type and concentration. It is not possible to give general rules for either of these variables and they have to be optimized experimentally. Important for the solubilization is the detergent-to-protein ratio. At low ratios (1:10) the membranes are lysed and large complexes of are formed containing protein, detergent, and membrane lipids. With progressively larger ratios smaller complexes are obtained. Finally, at ratios of 10:1 to 20:1 individual detergent-protein complexes are formed free of membrane lipids. To determine the optimal conditions it is important to vary both the detergent and the protein concentration.

[0028] Commonly used detergents and their critical micel concentrations (cmc) are listed in Table 3.

Table 3

Class of detergent	examples	cmc (mM)
Non-ionic	Triton X-100	0.30
	octylglucoside	25
Zwitterionic	CHAPS	6.5
	Zwittergent 3-12	3.6
Ionic	sodium deoxycholate	4 (at pH >8)

Isolation of inclusion bodies

[0029] In the cell there is competition between folding and aggregations. In many cases and in several host systems, recombinant proteins accumulate intracellularly in insoluble aggregates. The proteins in these so-called inclusion bodies are mostly inactive and denatured. In addition, dimers and multimers may be present. However, the expression of recombinant proteins in inclusion bodies can also be advantageous:

- the recombinant protein deposited in inclusion bodies can be 50% or more of the total cellular protein.
- the inclusion bodies often contain almost exclusively the overexpressed protein.
- in inclusion bodies the protein is protected from proteolytic degradation.
- expression in inclusion bodies will protect the cell against the toxicity of the recombinant protein.

The major problem is to recover biologically active and/or soluble protein in high yield. In order to accomplish this, the protein in the inclusion bodies must be solubilized and refolded *in vitro*. This procedure is carried out in three phases:

[0030] Inclusion bodies have a relatively high density and, therefore, can be pelleted by centrifugation. Cells are usually disrupted by high-pressure homogenization (optionally following a lysosyme treatment). It is important that cell lysis is complete, because intact cells sediment together with the inclusion bodies, thus contaminating the preparation. After centrifugation, the pellet is washed with buffer containing either low concentrations of chaotropic agents (e.g. 0.5 - 1 M guanidine-HCl or urea) or detergents (e.g. 1% Triton X-100 or 1 mg/ml sodium deoxycholate). This wash step is necessary to remove contaminants, especially proteins (proteases), that may have absorbed onto the hydrophobic inclusion bodies during processing.

Solubilization of the aggregated proteins in the inclusion bodies

[0031] The washed inclusion bodies are resuspended and incubated in buffer containing a strong denaturant and a reducing agent (usually 20 mM DTT or β -mercaptoethanol). The addition of a reducing agent keeps all cysteines in the reduced state and cleaves disulfide bonds formed during the preparation. Incubation temperatures above 30°C

are typically used to facilitate the solubilization process.

[0032] Optimal conditions for solubilization are protein specific and have to be determined for each protein.

[0033] After solubilization the solution can be centrifuged or filtered to remove remaining aggregates which could act as nuclei to trigger aggregation during refolding.

5 [0034] According to a specific embodiment of the invention reported here after in the Examples, this step is indicated as the Step 1 of the specific process described.

Refolding of the solubilized proteins

10 [0035] Refolding of the solubilized proteins is initiated by the removal of the denaturant. The efficiency of refolding depends on the competition between correct folding and aggregation. To slow down the aggregation process refolding is usually carried out at low protein concentrations, in the range of 10-100 mg/ml. Furthermore, refolding conditions must be optimized for each individual protein. Important variables are the buffer composition (pH, ionic strength) the temperature and the additives (often in combination).

15 [0036] If proteins contain disulfide bonds, the refolding buffer has to be supplemented with a redox system. The addition of a mixture of reduced and oxidized forms (1-3 mM reduced thiol and a 5:1 to 1:1 ratio of reduced to oxidized thiol) of low molecular weight thiol reagent usually provides the appropriate redox potential to allow formation and reshuffling of disulfide bonds. The most commonly used redox shuffling reagents are reduced and oxidized glutathione, but also cysteine and cysteamine are used.

20 [0037] For certain protein, probably due to low solubility of folding intermediates, this procedure is not very effective. Alternatively, the protein is completely oxidized in the presence of a large excess of oxidized glutathione, followed by dilution in refolding buffer containing catalytic amounts of reduced glutathione.

[0038] Different methods for the refolding of proteins have been described. The most used method is the removal of the solubilizing agent by dialysis. During dialysis the concentration of the solubilizing agent decreases slowly which 25 allows the protein to refold optimally. The ratio of the volumes of the sample and the dialysis buffer should be as such that at the equilibrium concentration of the solubilizing agent the protein has completely refolded.

[0039] The concentration of the solubilizing agent is decreased by dilution allowing the protein to refold. Usually the dilution is carried out slowly by step-wise addition of buffer or by continuous addition using a pump.

[0040] During dialysis and slow dilution the protein is exposed for an extended period of time to an intermediate concentration of the solubilizing agent (2-4 M urea or guanidine-HCl) where it is not yet folded but no longer denatured and thus extremely prone to aggregation. This could be prevented by the rapid dilution of the solubilized protein solution into the refolding buffer. Aggregation during this process can be limited by adding mild solubilizing agents to the refolding buffer, such as non-detergent sulfobetaines.

30 [0041] In order to keep the concentration of the unfolded protein low, thus limiting aggregation, aliquots of denatured protein are added at defined time points to the refolding buffer. The time intervals between two pulses have to be optimized for each individual protein. The process is stopped when the concentration of denaturant reaches a critical level with respect to refolding of the specific protein.

[0042] According to a specific embodiment of the invention reported here after in the Examples, this step is indicated as step 3 of the specific process described.

40 [0043] The solubilizing agent is removed using a chromatographic step. The application of different chromatography methods have been described:

- size exclusion chromatography (e.g. gel filtration on a Superdex 75 column)
- ion exchange chromatography
- 45 • affinity chromatography (e.g. IMAC using Chelating Sepharose or Ni-NTA agarose)

[0044] The denaturant is removed while the protein is slow migrating through the column or bound to the matrix. This usually gives a high yield of active protein even at protein concentrations in the mg/ml range.

50 EXAMPLES

[0045] The invention will now be described in detail with respect to the purification of a Chemokine expressed in E. Coli. Chemokines constitute a family of small pro-inflammatory cytokines with leukocyte chemotactic and activating properties. Depending on the position of the first conserved cysteines, the chemokine family can be divided in C-C, C-X-C and C-X₃-C chemokines (Baggiolini M. et al., Adv Immunol. 1994, 55:97-179; Baggiolini M. et al., Annu Rev Immunol. 1997, 15:675-705; Taub D. et al., Cytokine Growth Factor Rev. 1996, 7(4):355-76).

55 [0046] In particular, the detailed description that follows reports the purification of a triple mutant of human RANTES. This protein in the mature form (that is, after cleavage of the leader sequence MKKKWPR), has the sequence of SEQ

EP 1 439 191 A1

ID NO: 1 and will be referred to in the further part of this specification as RANTES triple mutant. The overall process is summarized by the flow chart presented in Figure 1, starting from the step of "solubilization of the aggregated proteins in the inclusion bodies". The two previous classical steps of preparation of cell lysates and isolation of inclusion bodies are carried out using the methods known in the art and/or as described before.

Step 1 - Solubilization of the Aggregated Proteins the Inclusion Bodies

[0047] 60-90 gr of inclusion bodies containing RANTES triple mutant are thawed. The pellet is dissolved with 450-800 ml of solubilization buffer (6M Guanidinium Chloride, 0.1 M Tris/HCl, 2mM DTT, pH 7.5+/-0.1) by means of the homogeniser Polytron until big particles are not visible anymore (about 5 minutes). Once homogenised, the solution is brought at a temperature of 60 ± 1 °C for 30'.

[0048] Then the solution is allowed to reach room temperature and filtered with a 1.2 µm membrane.

[0049] The material is usually processed as soon as it is prepared, in alternative it can be stored at -80°C.

Step 2 - Reverse Phase Chromatography On Source 30 RPC

Buffers

Pre-equilibration buffer: 0.1 M TRIS/HCl pH 7.5±0.1, 5% (v/v) Acetonitrile.

[0050] 12.1 g. of Tris/hydroxymethylaminomethane is added to 800-900 ml of purified water. The pH is adjusted to 7.5 with HCl. 50ml of Acetonitrile is added and the solution is brought to 1 liter with purified water.

Equilibration buffer: 0.1 M TRIS/HCl pH 7.5±0.1, 6M Guanidine and 2 mM DTT, 5% (v/v) Acetonitrile.

[0051] 6M guanidine and 2mM DTT are added to 600 ml of pre-equilibration buffer under stirring. After dissolution is completed, the solution is brought to 1 liter with pre-equilibration buffer.

Washing buffer: 5% (v/v) Acetonitrile, 0.1 % TFA in water

[0052] 50 ml of acetonitrile and a vial of TFA are added to 900 ml of purified water and the solution is brought to 1 liter with purified water.

Elution buffer: 35% (v/v) Acetonitrile, 0.1 % TFA in water

[0053] 350 ml of acetonitrile and a vial of TFA are added to 500-600 ml of purified water and the solution is brought to 1 liter with purified water.

Regeneration buffer: NaOH 0.5 M, 60% n-propanol

[0054] 20 gr of NaOH is added to 600 ml of n -propanol and the solution is brought up to 1 liter with purified water.

Storage buffer: NaOH 10 mM

[0055] 0.4 gr of NaOH is added to 1000 ml of purified water

Procedure

[0056] A column packed with Source 30RPC resin is flushed with at least 1 BV of NaOH 0.5M and subsequently with purified water until the pH reaches neutral values. Then the column is equilibrated with 3 BV of pre-equilibration buffer followed by 2-3BV of the equilibration buffer. The pH is checked and washing is continued if the parameters of the column's effluent are out of target values: pH 7.5 ± 0.1 .

[0057] The solubilized material prepared as described above is then loaded on the resin in the range of 5-18 mg of RANTES triple mutant/ ml of resin. When sample loading is completed, the column is flushed with 3-4 BV of equilibration buffer 0.1 M Tris/HCl + 5% Acetonitrile, pH 7.5 ± 0.1 . The flow trough and the wash are collected together.

[0058] Then the column is flushed with the washing buffer: 5% Acetonitrile + 0.1% TFA in water until the pH of the column effluent is acidic (control with the pH indicator).

[0059] Elution is started with the elution buffer: 35% Acetonitrile + 0.1% TFA in water. The fractions are collected as

EP 1 439 191 A1

soon as the OD signal starts to rise. About 3 -4 BV are collected as elution fraction.

[0060] The protein content of the fraction is calculated by the OD 280 nm analysis ($\epsilon=2.1$). The fraction is stored at $+5\pm 3^{\circ}\text{C}$ waiting for the next step of refolding.

[0061] After elution is completed flush the column with at least 4 BV of regeneration buffer followed by 3BV of purified water.

[0062] The column is flushed with at least 3 BV of NaOH 0.5M, and then the column is rinsed with purified water.

[0063] Afterwards the column is flushed with at least 3 BV of storage buffer and stored at room temperature until the next cycle.

[0064] The step here described is very efficient in the removal of cell impurities and gives a solution containing semipurified RANTES triple mutant ready to be submitted to the refolding step.

Step 3 - Refolding

[0065] In this step the unfolded RANTES triple mutant, eluted from the RPC column, is refolded by means of dilution in a suitable buffer and the use of a redox system.

Refolding buffer: 0.1 M Tris/HCl pH 7.5, 0.2 mM glutathione reduced form, 0.02 mM Glutathione oxidized form; cond. 6.0 ± 1 mS/cm

[0066] 12.1 gr of Tris, 61.5 mg of glutathione reduced and 12.2 mg of glutathione oxidized are added to 900 ml of purified water under stirring, bring to pH 7.5 ± 0.1 with concentrated HCl and make up to 1 liter. The solution is then stored at refrigerated temperature $+5^{\circ}\text{C} \pm 3^{\circ}\text{C}$ for not more than two days.

Procedure

[0067] This step is performed at refrigerated temperature ($+5\pm 3^{\circ}\text{C}$). The RPC elution fraction is diluted by adding it overnight drop by drop under mild stirring into a suitable volume of refolding buffer in order to reach a theoretical final concentration (calculated by OD) of 0.2-0.4 mg/ml.

[0068] Step 2 described before is very efficient in the removal of cell impurities giving a solution containing semipurified RANTES triple mutant ready to be submitted to the refolding step.

Step 4 - IEC on SP Sepharose FF

Buffers And Solutions

Equilibration buffer: 0.1 M Tris/HCl pH 7.5 ± 0.1 , conductivity 6.0 ± 1 mS/cm

[0069] 12.1 g of Tris is added to 900 ml of purified water under stirring, pH is brought to 7.5 ± 0.1 with 37% HCl and the solution is brought up to 1 liter. Conductivity is checked. The solution is stored at room temperature for not more than two days.

Elution buffer: 0.1 M Tris/HCl pH 7.5 ± 0.1 , 0.6 M NaCl, conductivity 57.0 ± 2.0 mS/cm

[0070] 12.1 g of Tris and 35 g of NaCl are added to 900 ml of purified water under stirring, the pH is brought to 7.5 ± 0.1 with HCl concentrated and the solution is brought up to 1 litre. Conductivity is checked. This solution is stored at room temperature $+20^{\circ}\text{C} \pm 5^{\circ}\text{C}$ and used within two days.

Regeneration buffer: 1.5 M NaCl

[0071] 87.6 g of NaCl is dissolved in 900 ml of purified water under stirring and the solution is brought up to 1 litre. This solution is stored at room temperature $+20^{\circ}\text{C} \pm 5^{\circ}\text{C}$ and used within two days.

Sanitizing solution: 0.5M NaOH

[0072] 20 g of NaOH is dissolved in 900 ml of purified water under stirring and the solution is brought up to 1 litre.

EP 1 439 191 A1

Storage solution: 0.01M NaOH

[0073] 0.4 g of NaOH is dissolved in 900 ml of purified water under stirring and the solution is brought up to 1 liter. This solution is stored at room temperature $+20\pm5^{\circ}\text{C}$ and used within two days.

Procedure

[0074] The column packed with SP Sepharose F resin is flushed with 3BV of NaOH 0.5 M followed by purified water until the pH goes down to neutral values (pH Indicator paper). Then the column is flushed again with 5 or more BV of equilibration buffer. pH and conductivity are checked and the wash is continued, if the parameters of the column's effluent are out of target values: pH 7.5 ± 0.1 , conductivity 6 ± 1 mS/cm.

[0075] Then the filtered post refolding material is loaded and the unbound fraction is collected.

[0076] When sample loading is completed, the column is flushed with 1 -2 BV of equilibration buffer and collection is continued. At this point usually the UV signal is on the baseline.

[0077] Then the column is flushed with the elution buffer. The elution fraction is started to be collected when the UV signal starts to rise, usually after the first 0.8 BV. 3BV of eluate is collected and this fraction is stored at $2-8^{\circ}\text{C}$ until the next step of cleavage.

[0078] After elution is completed, the column is flushed with at least 4 BV of regeneration buffer.

[0079] The column is flushed with at least 3 BV of NaOH 0.5M and then rinsed with 3 BV of purified water. The column is flushed with at least 3 BV of storage solution and stored until the next cycle.

Step 5 - Cleavage

[0080] The cleavage step is used to remove the -RANTES triple mutant leader sequence (MKKKWPR) from the correct sequence.

Procedure

[0081] The cleavage is performed at a temperature of $37\pm1^{\circ}\text{C}$ in a thermostatic oven.

[0082] The quantity of trypsin to be added is calculated on a ratio of 1:1 0000 in respect of the RANTES triple mutant content calculated by RP-HPLC analysis.

[0083] The trypsin is added as soon as the material to be cleaved has reached the target temperature.

[0084] The ongoing of the cleavage is followed by the RP -HPLC analysis that is able to separate the not cleaved peak from the cleaved one. Usually the cleavage is completed after about 3 hours. Just after the cleavage is completed, in order to stop the reaction, the pH of the solution is brought to 3.2 with concentrated phosphoric acid.

Step 6 - IEC on SP Sepharose HP

Buffers And Solutions

Buffer A -Equilibration-: 50 mM Ammonium Acetate pH 3.2 ± 0.2 , conductivity 0.4 ± 0.1 mS/cm

[0085] 2.8 ml of acetic acid is added to 900 ml of purified water under stirring, the pH is brought to 3.2 ± 0.2 with ammonia 25% and the solution is brought up to 1 litre.

Conductivity is checked. The solution is stored at room temperature $20^{\circ}\text{C}\pm5^{\circ}\text{C}$ for not more than two days.

Buffer B -Elution-: 50 mM Ammonium Acetate pH 3.2 ± 0.2 , 1 M NaCl conductivity 83 ± 2 mS/cm

[0086] 2.8 ml of acetic acid and 58.4 g of NaCl are added to 900 ml of purified water under stirring, the pH is brought to 3.2 ± 0.2 with ammonia 25% and the solution is brought up to 1 litre. pH and conductivity are checked.

[0087] The solution is stored at room temperature $20^{\circ}\text{C}\pm5^{\circ}\text{C}$ for not more than two days.

Sanitizing solution: 0.5M NaOH

[0088] 20 g of NaOH is dissolved in 900 ml of purified water under stirring and the solution is brought up to 1 litre.

EP 1 439 191 A1

Storage solution: 0.01 M NaOH

[0089] 0.4 g of NaOH is dissolved in 900 ml of purified water under stirring and the solution is brought up to 1 litre. This solution is stored at room temperature $+20\pm 5^{\circ}\text{C}$ and used within five days.

Procedure

[0090] Post cleavage material is diluted 1:2 with purified water in order to obtain a conductivity of the solution around 30 mS/cm

[0091] This chromatographic step requires the use of a system able to perform gradients between the equilibration and elution buffers.

[0092] The column is flushed with at least 3BV of NaOH 0.5 M and rinsed with purified water until the pH is neutral. Then the column is flushed with column 4-5 or more BV of equilibration buffer. pH and conductivity are checked and washing is continued, if the parameters of the column's effluent are out of target values:

pH 3.2 ± 0.2 and conductivity 0.4 ± 0.1 mS/cm

[0093] The starting material prepared as above is loaded on the column. As soon as the loading is completed, the column is flushed with 2-3 BV of equilibration buffer, until the UV signal reaches the baseline and then the column is flushed with a buffer composed by 60% of equilibration buffer and 40% of elution buffer. About 5-6 BV of this washing fraction is collected.

[0094] Elution is performed with a gradient between 40% and 100% of buffer B in 20 BV. The elution is collected in fractions (1 min. each one) that will be pooled later on, following the chromatogram pattern that usually shows the peaks from the peptide sequence, the cleaved (and correct) one and the uncleaved molecule. Based on the chromatogram the fractions that form the peak of interest are collected in order to form the elution fraction.

[0095] Then the column is flushed with at least 3-4 BV of NaOH and then rinsed with 3-4 BV of purified water.

[0096] The column is flushed with at least 3 BV of storage solution and stored until the next cycle.

Step 7 - Bulk Ultrafiltration

Buffers And Solutions

Buffer: 50 mM Ammonium Acetate, pH 4.0 ± 0.1 conductivity 1.0 ± 0.2 mS/cm

[0097] 2.8 ml of Acetic Acid is added to 900 ml of purified water, under stirring. The solution is brought up to 1 litre. The pH is brought to 4.0 ± 0.1 with 25% ammonia and conductivity is checked. This solution is stored at room temperature $+20\pm 5^{\circ}\text{C}$ and used within 1 day.

Sanitizing solution: 0.5 M NaOH

[0098] 20 g of sodium hydroxide are dissolved in 900 ml of purified water under stirring. The solution is brought up to 1 litre and stored at room temperature.

Storage solution: 0.05 M NaOH

[0099] 2 g of sodium hydroxide is dissolved in 900 ml of purified water under stirring, the solution is brought up to 1 litre, stored at room temperature $+20\pm 5^{\circ}\text{C}$ and used within three days.

Procedure

[0100] The ultrafiltration system is assembled with a regenerated cellulose membrane at cut off 3KD (Millipore).

[0101] The ultrafilters are sanitized with 200 ml of NaOH 0.5 M by recycling the alkaline solution for not less than 30 minutes and then rinsed with purified water until permeate pH is below 7.5.

[0102] The SP HP pool elution solution is recirculated on the system at a constant pressure of slightly less than 0.5 bar and purified water is added, trying to maintain the volume constant. Conductivity is checked and washing is continued until it reaches values below 1 mS/cm.

[0103] Then the column is washed with the bulk buffer: 50 mM Ammonium Acetate pH 4.0 ± 0.1 , cond. 1.0 ± 0.2 mS/cm until permeate reaches the same values.

EP 1 439 191 A1

[0104] The retentate fraction is collected and the system is washed recirculating it with the same buffer until no foam is visible anymore.

[0105] The ultrafilters are washed and sanitized by recycling around 200 ml of 0.5 M NaOH for at least 30 minutes, then rinsed ultrafilters with purified water until the permeate pH is below 7.5. The ultrafilters are then stored in 0.05M NaOH at +20 ±5°C until the next cycle.

Analysis

[0106] The preliminary analysis for quantitation of RANTES triple mutant has been set up to measure the concentration of the molecule in the in process samples to monitor the yield of the process and to monitor both refolding and cleavage steps. Here follows the description of the analyses.

Equipment	Analytical HPLC system equipped with column oven		
Column	Symmetry C18, 3.5 µ, 4.6x 75mm		
Eluent A	0.1% aqueous TFA		
Eluent B	0.1% TFA in acetonitrile		
Eluent C (for column storage)	50% acetonitrile in water		
Method			
UV detection	214 nm		
Temperature	45°C		
Injection time	50 minutes		
Injection volume	10-100 µl		
Gradient	Time	%A	%B
	0	75	25
	25'	65	35
	25.1'	20	80
	30'	20	80
	*30.1'	75	25
	*38'	75	25
	*38.1'	75	25
Flow rate	1 ml/minute * step performed at 1.2 ml/min		
Reference sample	Purified uncleaved RANTES triple mutant at the concentration of 0.65 g/ml by OD (E=21)		

[0107] The overall recovery of the purification process estimated by RP-HPLC resulted to be over 30%, which is a very good result compared with the data in the literature; the purity by RP-HPLC is >90%, by SE-HPLC is > 95% and the HCP content, analysed with an ELISA commercial kit, is < 100 ppm.

[0108] A panel of such other analyses as SDS, IEF and Maldi TOF on the several batches of drug substances produced confirmed the consistency of the process both as purity and quality of the molecule.

[0109] Although the process has been here described with abundance of particulars with respect to the purification of a specific chemokine mutant, it is understood that many obvious variations are within the skill of an average expert in the field.

[0110] Such variations of the process parameters, as well as the application of the process to any protein that is able to withstand the process conditions, may be made without affecting the spirit and scope of the invention that is defined by the following claims.

SEQUENCE LISTING

5 <110> APPLIED RESEARCH SYSTEMS ARS HOLDING N.V.
 <120> PROCESS FOR THE PURIFICATION OF BACTERIALLY EXPRESSED PROTEINS
 <130> EP808Z
 10 <160> 2
 <170> PatentIn version 3.1
 <210> 1
 15 <211> 67
 <212> PRT
 <213> Escherichia coli
 <400> 1
 20 Ser Pro Tyr Ser Ser Asp Thr Thr Pro Cys Cys Phe Ala Tyr Ile Ala
 1 5 10 15
 25 Arg Pro Leu Pro Arg Ala His Ile Lys Glu Tyr Phe Tyr Thr Ser Gly
 20 25 30
 Lys Cys Ser Asn Pro Ala Val Val Phe Val Thr Ala Ala Asn Ala Gln
 30 35 40 45
 Val Cys Asn Pro Glu Lys Lys Trp Val Arg Glu Tyr Ile Asn Ser Leu
 50 55 60
 35 Glu Met Ser
 65
 40 <210> 2
 <211> 74
 <212> PRT
 <213> Escherichia coli
 45 <220>
 <221> MISC_FEATURE
 <222> (1)..(7)
 <223> Leader sequence
 50
 55

<400> 2

5 Met Lys Lys Lys Trp Pro Arg Ser Pro Tyr Ser Ser Asp Thr Thr Pro
1 5 10 15

10 Cys Cys Phe Ala Tyr Ile Ala Arg Pro Leu Pro Arg Ala His Ile Lys
20 25 30

15 Glu Tyr Phe Tyr Thr Ser Gly Lys Cys Ser Asn Pro Ala Val Val Phe
35 40 45

Val Thr Ala Ala Asn Ala Gln Val Cys Asn Pro Glu Lys Lys Trp Val
50 55 60

20 Arg Glu Tyr Ile Asn Ser Leu Glu Met Ser
65 70

25 2

Claims


- 30 1. In a process for the recovery of a protein expressed in prokaryotic host cells as inclusion bodies and its subsequent purification, the improvement consisting in interposing a Reverse Phase Chromatography step between the step of solubilization of the aggregated proteins in the inclusion bodies/denaturation and the renaturation/refolding step.
- 35 2. The process of claim 1, in which the prokaryotic cells are bacterial cells.
3. The process of claims 1 or 2, in which the bacterial cells are E. Coli cells.
- 40 4. An improved process according to any of claims 1 and 2, in which the protein extracted and purified is a chemokine.
5. The process of claim 3 wherein the chemokine is the chemokine mutant of SEQ ID NO: 1.

45

50

55

Figure 1

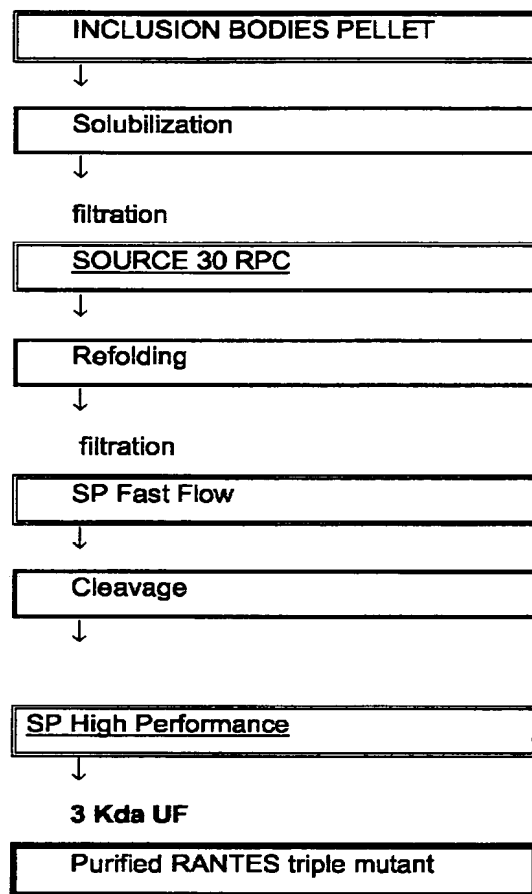



European Patent
Office

EUROPEAN SEARCH REPORT

Application Number
EP 03 10 0064

DOCUMENTS CONSIDERED TO BE RELEVANT			
Category	Citation of document with indication, where appropriate, of relevant passages	Relevant to claim	CLASSIFICATION OF THE APPLICATION (Int.Cl.7)
X	MUSACCHIO A ET AL.: "Recombinant Opc meningococcal protein, folded in vitro, elicits bactericidal antibodies after immunization" VACCINE, vol. 15, no. 6/7, 1996, pages 751-758, XP004064539 * page 752 *	1-3	C07K14/52 C12N15/19
A	WO 98 14467 A (HUMAN GENOME SCIENCES, INC.) 9 April 1998 (1998-04-09) * page 23 - page 27; example 1 *	4,5	
A	WO 02 28419 A (APPLIED RESEARCH SYSTEMS ARS HOLDING N.V.) 11 April 2002 (2002-04-11) * claims 7,14; example 1B *	4,5	
			TECHNICAL FIELDS SEARCHED (Int.Cl.7)
			C07K C12N
The present search report has been drawn up for all claims			
Place of search THE HAGUE		Date of completion of the search 15 April 2003	Examiner Schmidt, Harald
CATEGORY OF CITED DOCUMENTS X : particularly relevant if taken alone Y : particularly relevant if combined with another document of the same category A : technological background O : non-written disclosure P : intermediate document T : theory or principle underlying the invention E : earlier patent document, but published on, or after the filing date D : document cited in the application L : document cited for other reasons : member of the same patent family, corresponding document			

EPO FORM 1503 (03.02.99) (p.4/6)

BEST AVAILABLE COPY

**ANNEX TO THE EUROPEAN SEARCH REPORT
ON EUROPEAN PATENT APPLICATION NO.**

EP 03 10 0064

This annex lists the patent family members relating to the patent documents cited in the above-mentioned European search report. The members are as contained in the European Patent Office EDP file on. The European Patent Office is in no way liable for these particulars which are merely given for the purpose of information.

15-04-2003

Patent document cited in search report		Publication date		Patent family member(s)	Publication date
WO 9814467	A	09-04-1998	US	6001606 A	14-12-1999
			US	5912327 A	15-06-1999
			EP	0975657 A1	02-02-2000
			JP	2001503614 T	21-03-2001
			WO	9814467 A1	09-04-1998
			US	2002061551 A1	23-05-2002
			US	6451562 B1	17-09-2002
			US	2002076746 A1	20-06-2002
WO 0228419	A	11-04-2002	US	6495129 B1	17-12-2002
WO 0228419	A	11-04-2002	AU	1591902 A	15-04-2002
			WO	0228419 A2	11-04-2002

EPO FORM P0459

For more details about this annex : see Official Journal of the European Patent Office, No. 12/82

BEST AVAILABLE COPY